


PRODUCT LATEX FOLEY CATHETERS

Raw Material:	made of siliconized latex
Tips:	rounded off, firm and safe, bending-proof, assuring an easy introduction. Available also Tiemann Tip (curved), which is designed to deal with the prostatic curve in male patients.
Eyes:	the big smooth eyes, facing each other, allow to minimize obstructions and to maximize drainage.
Way:	2 or 3 ways
Length :	for adult 410 mm.
Dimensions (CH):	Paediatric: CH 6,8, 10. Adult: from CH 12 to CH 24.
Balloon:	3 or 5 ml for paediatric catheters, 5/15 ml for others catheters, it is fully inspected by filling it with water and then by deflating it so as to guarantee its quality and proper working.
Catheter shaft:	the presence of silicone in the mix allows to prevent torsion and to minimize encrustments as well as traumas.
Y-junction:	thick and reinforced, it prevents contractions and torsions.
Funnel:	provided with special valve for Luer and Luer-Lock syringes, it guarantees a reliable inflating and deflating of the balloon while preventing possible leaks.
Packaging:	each single catheter is sealed into an inner sheath and then packed into a peelable bag in order facilitate its proper aseptic moving and to offer the highest guarantee of its sterility.
Sterilization:	made by ETO and CO2 blend, expiry date from the sterilization: 5 years.
Registered TradeMark:	Securdrain Vià!

DRAWING

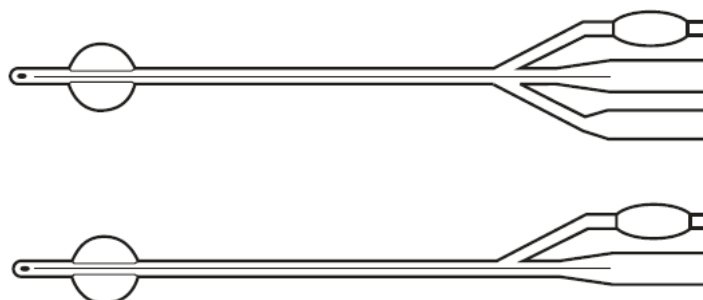

CODES AND PACKAGE

CODE	DESCRIPTION	PACKAGE
11601312061	Latex Foley Catheter, pediatric, 2 ways, 28 cm, ball 3 ml sterile, CH 06	10 pcs / box - 200 pcs/ carton
11601312082	Latex Foley Catheter, pediatric, 2 ways, 28 cm, ball 5 ml sterile, CH 08	10 pcs / box - 200 pcs/ carton
11601312103	Latex Foley Catheter, pediatric, 2 ways, 28 cm, ball 5 ml sterile, CH 10	10 pcs / box - 200 pcs/ carton
11601322121	Latex Foley Catheter, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 12	10 pcs / box - 200 pcs/ carton
11601322141	Latex Foley Catheter, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 14	10 pcs / box - 200 pcs/ carton
11601322161	Latex Foley Catheter, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 16	10 pcs / box - 200 pcs/ carton
11601322181	Latex Foley Catheter, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 18	10 pcs / box - 200 pcs/ carton
11601322201	Latex Foley Catheter, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 20	10 pcs / box - 200 pcs/ carton
11601322221	Latex Foley Catheter, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 22	10 pcs / box - 200 pcs/ carton
11601322241	Latex Foley Catheter, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 24	10 pcs / box - 200 pcs/ carton

11601363161	Latex Foley Catheter, 3 ways, 41 cm, balloon 5 -15 ml sterile, CH 16	10 pcs / box - 200 pcs/ carton
11601353161		
11601363181	Latex Foley Catheter, 3 ways, 41 cm, balloon 5 -15 ml sterile, CH 18	10 pcs / box - 200 pcs/ carton
11601353181		
11601363201	Latex Foley Catheter, 3 ways, 41 cm, balloon 5 -15 ml sterile, CH 20	10 pcs / box - 200 pcs/ carton
11601353201		
11601363221	Latex Foley Catheter, 3 ways, 41 cm, balloon 5 -15 ml sterile, CH 22	10 pcs / box - 200 pcs/ carton
11601353221		
11601363241	Latex Foley Catheter, 3 ways, 41 cm, balloon 5 -15 ml sterile, CH 24	10 pcs / box - 200 pcs/ carton
11601353241		

11601342121	Latex Foley Catheter Tiemann Tip, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 12	10 pcs / box - 200 pcs/ carton
11601342141	Latex Foley Catheter Tiemann Tip, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 14	10 pcs / box - 200 pcs/ carton
11601342161	Latex Foley Catheter Tiemann Tip, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 16	10 pcs / box - 200 pcs/ carton
11601342181	Latex Foley Catheter Tiemann Tip, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 18	10 pcs / box - 200 pcs/ carton
11601342201	Latex Foley Catheter Tiemann Tip, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 20	10 pcs / box - 200 pcs/ carton
11601342221	Latex Foley Catheter Tiemann Tip, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 22	10 pcs / box - 200 pcs/ carton
11601342241	Latex Foley Catheter Tiemann Tip, 2 ways, 41 cm, balloon 5 -15 ml sterile, CH 24	10 pcs / box - 200 pcs/ carton